



Many dental treatments and procedures start with the taking of impressions. High quality impressions give the dentist confidence in an exceptional end result. Cavex offers a tried and tested line of silicones and auxiliary materials developed according to the requirements of the dental practice.

*Your impression is our concern*

## Cavex SiliconA bite, perfect the first time – every time

A good bite registration should be a pleasure to work with, enabling the registration to be perfect the first time – every time. This can only be possible when the material is easy to handle, easy to trim and absolutely accurate.

Cavex SiliconA bite “at least” meets these criteria by having a short working and setting time, both 30 seconds, and having a hardness measured at a Shore D of 40. This hardness is achieved by having a highly filled A-silicone material, enabling the bite registration to be 100% accurate.

The enclosed Bite Registration Tips ensure a consistent, even layer of material completely covering the surface of the elements, enabling a perfect bite registration.






**Cavex Holland BV**

P.O. Box 852  
2003 RW Haarlem  
The Netherlands

T +31 23 530 77 00  
F +31 23 535 64 82  
dental@cavex.nl  
www.cavex.nl

## Every detail registered

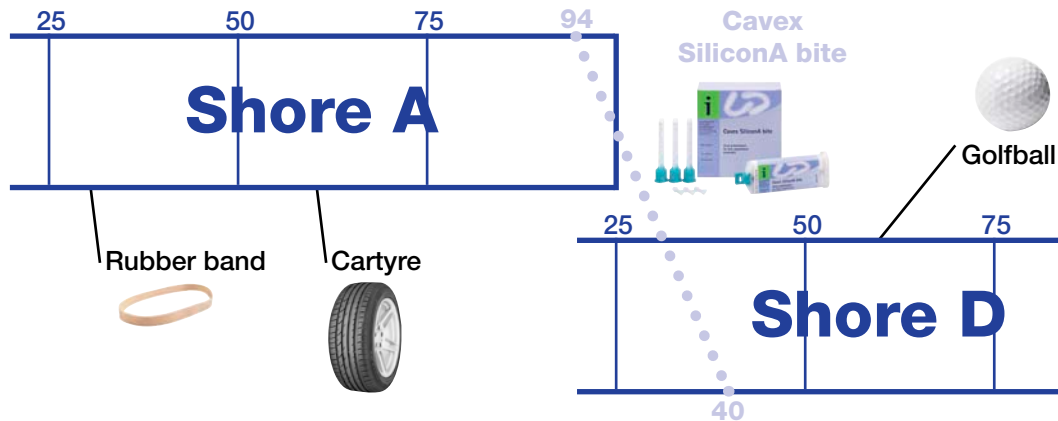


- **Superior hardness** Cavex SiliconA bite's hardness of Shore D 40, equal to a Shore A of >94, ensures a perfectly accurate registration of the occlusion which stays exactly the same through time. The material is easy to trim with a scalpel or grinder, which makes it a joy to work with for the Dental Technician.
- **Fast** The working and setting time ensure a short chair-time since both take just 30 seconds. Therefore the whole procedure takes only 60 seconds, reducing the patients discomfort to a minimum. The fast 30 second setting time minimizes the chances of discrepancy.
- **Highly thixotropic** An important stage of the procedure is applying the material on the teeth. The Bite Registration Tips enable a complete, even layer of Cavex SiliconA bite to be applied on the teeth. The high thixotropic nature of the material makes sure it stays in the right place. When dispensing the material, it flows smoothly out of the cartridge. Once applied to the teeth it stops flowing, ruling out the danger of uncovered teeth or an uneven layer.
- **Scannable** Cavex SiliconA bite is such a high filled material, that it is perfectly suitable to scan, without using powder.

For further information, please contact Cavex: +31 23 530 77 00 or [www.cavex.nl](http://www.cavex.nl)

- **Unique dispensing tips**



bite registration tip

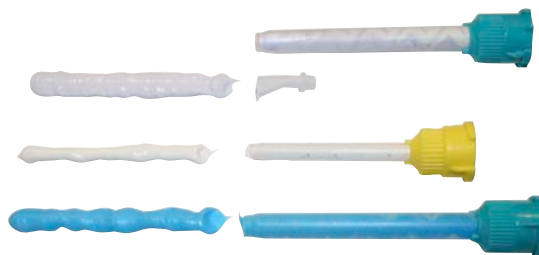

different dispensing tips

### Assortment:

**A554** Cavex SiliconA bite

2 x 50ml cartridges